

javelin

The Javelin uses the best features of composite technology to offer comfortable, lively and secure gait action. The foot is streamlined to the basic elements that govern walking and the result: good flexibility and support, excellent terrain compliance and a super lightweight walking experience.

- The Javelin offers a wide size range from 22cm to 30cm, with a choice of light or dark tone foot shells
- The slim-line continuous shin and toe spring flexes to absorb shock and reduce strain on the back and muscles
- Designed for the lifestyles of K2 and K3 ambulators
- 365mm max build height 200mm min build height
- Minimal design offers easy solution to complementary cosmetic finishing
- 36 month warranty



Part number

example:

JAV 22L 3

Size/side Spring set

If a dark tone foot shell is required please add suffix "D" to the part number.
The top housing is not supplied with the foot. Please order separately.
The fitting clamp can be used for alignment and set up, to ensure correct length before cutting the shin.

Spring Selection

Activity	Amputee Weight (kg)									
	100-115 44-52	116-130 53-59	131-150 60-68	151-170 69-77	171-195 78-88	196-220 89-100	221-255 101-116	256-285 117-130	286-325 131-147	326-356 148-166
K2	1	1	2	3	4	5	6	7	8	9
K3	1	2	3	4	5	6	7	8	9	

Spring set

A walk in the park

can be

a

STROLL

or an

ADVENTURE.


For great

off-road

RESPONSE

try

JAVELIN.

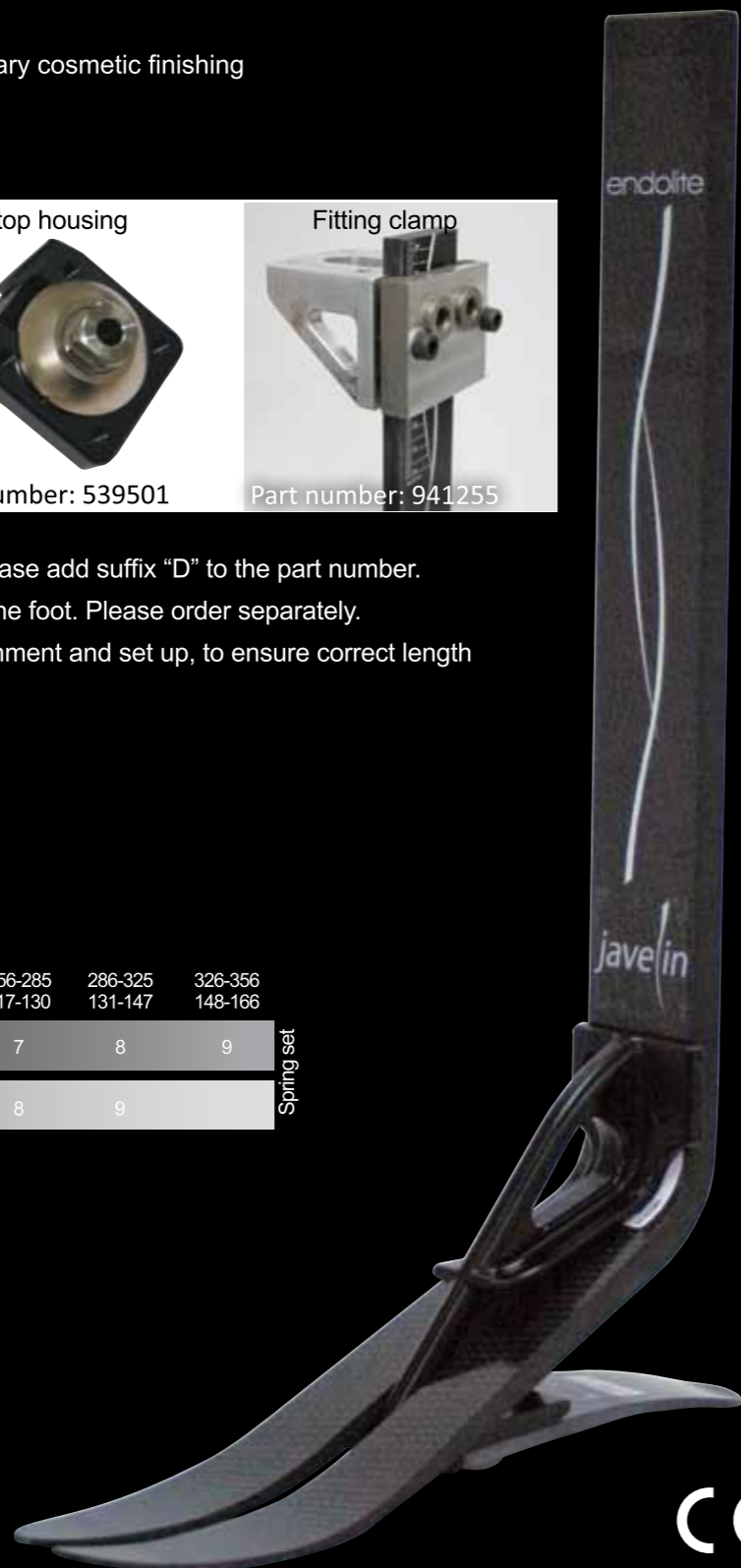

javelin

The inscrutable Sphinx famously set a riddle for all travellers into Thebes, "What walks on four feet in the morning, two in the afternoon and three at night?" This was finally answered by Oedipus, thus saving himself from being eaten by the Sphinx! "Man; as an infant crawls on all fours, as an adult, walks on two legs and, in old age, relies on a walking stick."

Today, thankfully, technology like the Javelin foot prolongs stable and independent walking, over all terrain, thus making encounters with a Sphinx a mere walk in the park!



The e-carbon heel gives supple control and excellent compliance



The active pylon flexes in a fluid and integrated kinematic response



Javelin's lively toe enhances the naturalness of individual gait



Slim good looks are not only the foot style but the cosmetic finish