
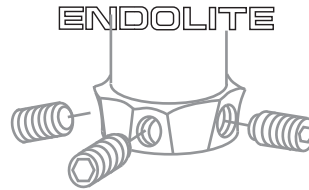


# ENDOLITE




**Blatchford Products Ltd**

Lister Road, Basingstoke, Hampshire, RG22 4AH, United Kingdom  
Telephone: +44 (0)1256 316600 Facsimile: +44 (0)1256 479705  
Email: sales@blatchford.co.uk www.blatchford.co.uk

## Super SACH Foot P529101 - P529114

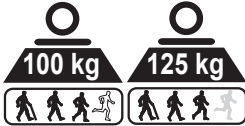

Adjuster screw  
maximum heel  
stiffness achieved  
with screw in fully  
engaged position  
(Do not over tighten)  
For minimum  
stiffness  
release screw  
by 4 turns

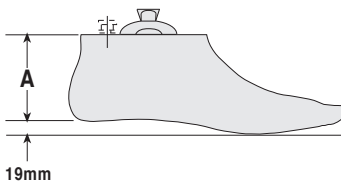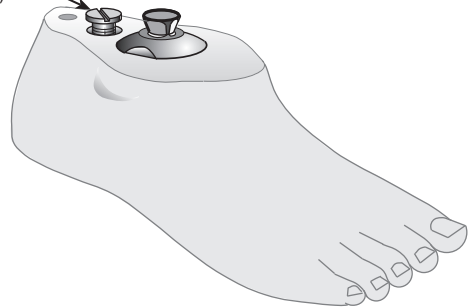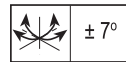

Size	A
22cm	} 70mm
23cm	
24cm	
25cm	} 75mm
26cm	
27cm	} 80mm
28cm	



938086/06-0705