

**Blatchford
Products Ltd**

Lister Road, Basingstoke, Hampshire,
RG22 4AH, United Kingdom
Telephone: +44 (0)1256 316600
Facsimile: +44 (0)1256 479705
Email: sales@blatchford.co.uk
www.endolite.co.uk




**Multiflex
Special
Edition
519525-532**



endolite
get busy living



***Before bonding,
lightly abraid and
clean mating
surfaces with
suitable solvent
cleaner.**

Contact
adhesive Evo-
stick 528
(926213),
Barge Cement
or equivalent*

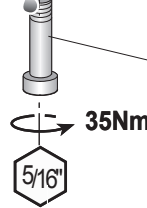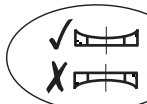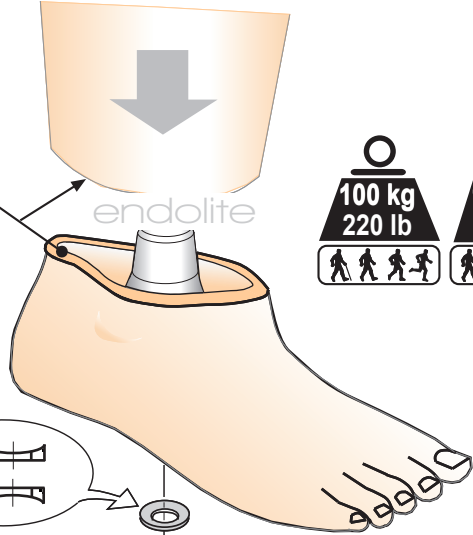

400022
A new bolt must be used when replacing the foot.

*Verify
engagement of
serrations
prior to
tightening foot
bolt*

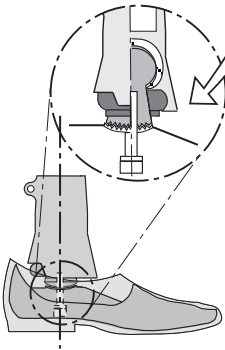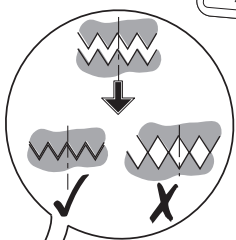

Size	A
25cm	} 75mm
26cm	
27cm	
28cm	} 80mm
28cm	

