

New Endolite Silicone Liners



Follow Blatchford: | blatchford.co.uk
@BlatchfordGrp | 0114 263 7900
/BlatchfordGrp | sales@blatchford.co.uk

Endolite is part of the Blatchford Group. Registered address, Chas A Blatchford & Sons Ltd, Unit D, Antura, Kingsland Business Park, Basingstoke, RG24 8PZ, United Kingdom. t: +44 (0) 1256 316600 f: +44 (0) 1256 316710



936433 lss2 05/16. Information correct at time of print.



Endolite Interface Technology,
know how that embodies comfort



Endolite Interface Philosophy

The increasing sophistication of prosthetic limb systems means that the interface is now even more critical to a successful outcome. Comfort, control and feedback parameters extend to the interactions between the body and the prosthesis. Endolite is proud to be developing technology in this area to enhance appropriate interface interactions and support seamless and comfortable prosthetic wearability for all users.

The Endolite Comfort Liner is gentle on fragile skin and protects residual limbs with areas susceptible to shear force abrasion and pressure peaks. The fabric cover with integrated matrix contours perfectly to the limb and zonal stretch areas allow comfortable knee flexion but resist distal stretch which would result in the prosthesis slipping away from the body.


The Cushion and Locking options allow the ideal choice of proximal or distal pin lock suspensions.

KEY FEATURES:



Stretch Zones
Proximal bi-directional stretch zone improves comfort and contouring and reduces shear force on the knee during flexion

Tendresse™ Silicone Finish
Reduced shear stress on skin because the Tendresse™ silicone finish means liners have a low coefficient of friction, in tests this measured 78% less than standard silicone

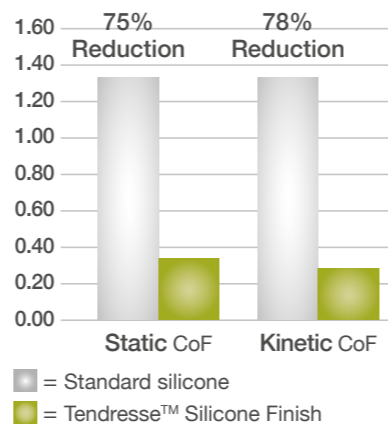

Extra Length Liner
At 485mm the Endolite Silicone Liners allow additional consideration for positioning the trim line

Integrated Matrix
Technical knitted outer cover for the Locking Liner has integrated distal matrix with uni-directional stretch to resist pistoning. The seamless knitting technology allows a smoother transition between the bi-directional and uni-directional stretch zones

Flexi-umbrella
The Locking Liner features a super flexi-umbrella which inverts easily to facilitate optimal donning



LINER SELECTION INFORMATION

Amputation level: Trans-tibial

Activity level: Medium to high


Sizes available: 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 40

Part Number example:

Endolite Comfort Cushion CTTCPXXL

Endolite Comfort Locking CCTLPPXXL

C	TT	CP/LP	XX	L	*Note: 23.5 = 23 26.5 = 26
Comfort	Trans-tibial	Cushion Parallel/ Locking Parallel	Size 22cm to 40cm*	Light Colour	

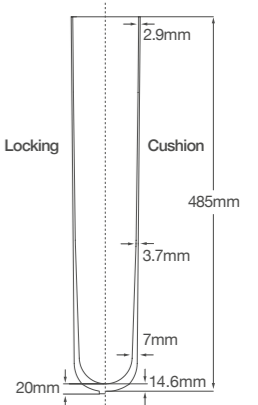

MEASUREMENT

Circumference taken 4cm from distal end of residual limb. Choose a smaller liner if the measurement is between sizes to ensure a snug fit.

