

Patented quality Made in Germany

# SoftSkin Softskin Air®

Get more information on  
[www.uniprox.de](http://www.uniprox.de)



The Liner Standard for above  
and below knee amputees.

unique prosthetic solutions  
A company of the Bauerfeind Group

 **uniprox**®

Patented quality Made in Germany

# Softskin Air® *No Sweat.*

Get more information on  
[www.uniprox.de](http://www.uniprox.de)



**Nick Weihs**  
Vice World Champion  
athletics

unique prosthetic solutions

 **uniprox®**



# Features SoftSkin and SoftSkin Air® Liners

- Use of medical silicone
- Available in 19 sizes for customization
- Unique functional and progressive outside texture
- Reduction of elongation
- Machine-washable (30°C)
- Breathable (SoftSkin Air®)
- Body moisture draining (SoftSkin Air®)
- Available as Trans Tibial (BK) and Trans Femoral (AK) version
- Available with or without distal connector

## SoftSkin

*New generation Silicone-Liner*

Silicone liners have proven to be accepted as the world's most successful suspension system. The quality and function of the Uniprox SoftSkin liners are known for more than ten years.

Thanks to the unique progressive outside textile, in combination with the medical silicone, a huge number of satisfied amputees uses the SoftSkin liners. The very durable and progressive textile provides constant compression over the entire length of the residual limb. The special mesh structure at the distal end, minimizes elongation. What makes it unnecessary for the SoftSkin liners to have additional reinforcements (matrix) inside the silicone.

## SoftSkin Air®

*The Breathable Liner – Tradition meets Innovation*

The SoftSkin Air® Liner is the world's first breathable Silicone Liner. Special micro pores in the silicone transport moisture (sweat) from the skin away through the liner wall. The results are a dryer skin and a lower temperature inside the liner. The additional effect is a better adhesion between the skin and the silicone to increase safety for the user.

Not only patients with sweating problems benefit from the features of the SoftSkin Air® liner, the liner allows entrapped air to escape through the micro pores what reduces risks of blisters.




# SoftSkin S30

The gentle liner



- Shore hardness: 30 / Thickness: 3 mm / Length: 370 mm
- Soft cushioning of the stump
- Good gripping effect of the silicone
- Comfortable and easy to handle
- For patients with mobility levels 1–3
- Easy donning due to special coating of the textile
- New textile cover provides good containment (no matrix needed)
- Stainless steel pin fitting with distal connector
- Inclusive liner stand

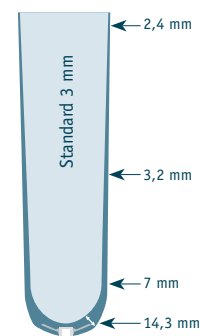

| Order No.                                | Deliverable sizes  | Description              |
|--|--|--------------------------|
| S30- <input type="text" value="Size"/>   | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | with distal connector    |
| S30CL- <input type="text" value="Size"/> | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | without distal connector |

Example: SoftSkin S30CL in Size 25 → Order No. S30CL-25


› The circumference is measured 4 cm above the distal end of the residual limb!

# SoftSkin Air® S30

The gentle breathable liner



- Shore hardness: 30 / Thickness: 3 mm / Length: 370 mm
- Soft cushioning of the stump
- Less sweating due to special micropores
- Long time wearing comfort
- Increased safety due to better friction on the skin
- Easy donning due to special coating of the textile
- New textile cover provides good containment (no matrix needed)
- Stainless steel pin fitting with distal connector
- Inclusive liner stand

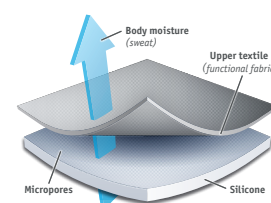

| Order No.                                     | Deliverable sizes  | Description              |
|---|--|--------------------------|
| S30- <input type="text" value="Size"/> -Air   | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | with distal connector    |
| S30CL- <input type="text" value="Size"/> -Air | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | without distal connector |

Example: SoftSkin Air® S30 in Size 36 → Order No. S30-36-Air


› The circumference is measured 4 cm above the distal end of the residual limb!

# SoftSkin S40

The comfortable liner



- Shore hardness: 40 / Thickness: 3 mm / Length: 370 mm
- Soft cushioning of the stump
- Good gripping effect of the silicone
- Comfortable and easy to handle
- For patients with mobility levels 1–3
- Easy donning due to special coating of the textile
- New textile cover provides good containment (no matrix needed)
- Stainless steel pin fitting with distal connector
- Inclusive liner stand

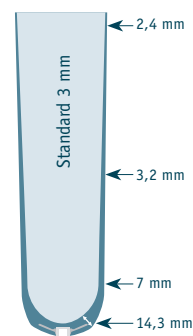

| Order No.                                | Deliverable sizes  | Description              |
|--|--|--------------------------|
| S40- <input type="text" value="Size"/>   | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | with distal connector    |
| S40CL- <input type="text" value="Size"/> | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | without distal connector |

Example: SoftSkin S40CL in Size 25 → Order No. S40CL-25


› The circumference is measured 4 cm above the distal end of the residual limb!

# SoftSkin Air® S40

The comfortable breathable liner



- Shore hardness: 40 / Thickness: 3 mm / Length: 370 mm
- Soft cushioning of the stump
- Less sweating due to special micropores
- Long time wearing comfort
- Increased safety due to better friction on the skin
- Easy donning due to special coating of the textile
- New textile cover provides good containment (no matrix needed)
- Stainless steel pin fitting with distal connector
- Inclusive liner stand

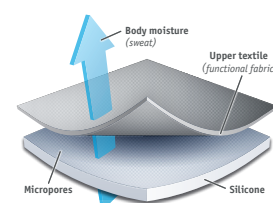

| Order No.                                     | Deliverable sizes  | Description              |
|---|--|--------------------------|
| S40- <input type="text" value="Size"/> -Air   | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | with distal connector    |
| S40CL- <input type="text" value="Size"/> -Air | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | without distal connector |

Example: SoftSkin Air® S40 in Size 36 → Order No. S40-36-Air


› The circumference is measured 4 cm above the distal end of the residual limb!

# SoftSkin S50

The active liner



- Shore hardness: 50 / Thickness: 3 mm / Length: 370 mm
- Hard silicone stabilizes the soft tissue
- Combination of firm silicone and textile provides optimal control longer
- For patients with mobility levels 2–4
- Easy donning due to special coating of the textile
- New textile cover provides good containment (no matrix needed)
- Stainless steel pin fitting with distal connector
- Inclusive liner stand

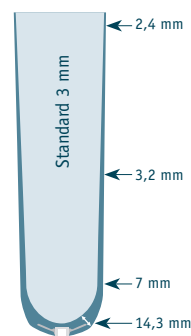

| Order No.                                | Deliverable sizes  | Description              |
|--|--|--------------------------|
| S50- <input type="text" value="Size"/>   | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | with distal connector    |
| S50CL- <input type="text" value="Size"/> | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | without distal connector |

Example: SoftSkin S50CL in Size 25 → Order No. S50CL-25


› The circumference is measured 4 cm above the distal end of the residual limb!

# SoftSkin Air® S50

The active breathable liner



- Shore hardness: 50 / Thickness: 3 mm / Length: 370 mm
- Hard silicone stabilizes the soft tissue
- Less sweating due to special micropores
- Long time wearing comfort
- Increased safety due to better friction on the skin
- Combination of firm silicone and textile = optimal control longer
- Easy donning due to special coating of the textile
- New textile cover provides good containment (no matrix needed)
- Stainless steel pin fitting with distal connector
- Inclusive liner stand

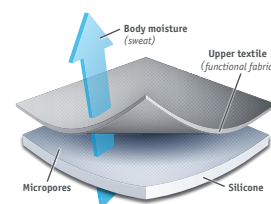

| Order No.                                     | Deliverable sizes  | Description              |
|---|--|--------------------------|
| S50- <input type="text" value="Size"/> -Air   | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | with distal connector    |
| S50CL- <input type="text" value="Size"/> -Air | 16 / 18 / 20 / 21 / 22 / 23,5 / 25 / 26,5 / 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 | without distal connector |

Example: SoftSkin Air® S50 in Size 36 → Order No. S50-36-Air

› The circumference is measured 4 cm above the distal end of the residual limb!

# SoftSkin TF

The comfortable liner



- Shore hardness: 40 / Thickness: 3 mm / Length: 350 mm
- Hard silicone stabilizes the soft tissue
- Soft cushioning of the stump
- Has high transverse elasticity so flexibly adapts to residual limb shape
- The special combination of fabrics effectively minimizes longitudinal travel
- For patients with mobility levels 1–4
- Easy donning due to special coating of the textile
- Stainless steel pin fitting
- Inclusive liner stand



| Order No.                              | Deliverable sizes                                    | Description           |
|--|--|-----------------------|
| STF- <input type="text" value="Size"/> | 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 / 50 / 55 | with distal connector |

Example: SoftSkin TF in Size 45 → Order No. STF-45


› The circumference is measured 4 cm above the distal end of the residual limb!

# SoftSkin Air® TF

The comfortable breathable liner



- Shore hardness: 40 / Thickness: 3 mm / Length: 350 mm
- Soft cushioning of the stump
- Less sweating due to special micropores
- Long time wearing comfort
- Increased safety due to better friction on the skin
- Has high transverse elasticity so flexibly adapts to residual limb shape
- The special combination of fabrics effectively minimizes longitudinal travel
- Stainless steel pin fitting
- Inclusive liner stand

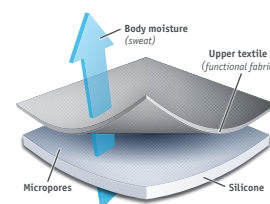

| Order No.                                   | Deliverable sizes                                    | Description           |
|---|--|-----------------------|
| STF- <input type="text" value="Size"/> -Air | 28 / 30 / 32 / 34 / 36 / 38 / 40 / 42 / 45 / 50 / 55 | with distal connector |

Example: SoftSkin Air® TF in Size 36 → Order No. STF-36-Air

› The circumference is measured 4 cm above the distal end of the residual limb!

Patented quality Made in Germany

# SoftSkin Softskin Air®

Get more information on  
[www.uniprox.de](http://www.uniprox.de)

**Uniprox GmbH & Co. KG**

Heinrich-Heine-Straße 4

07937 Zeulenroda

Germany

Phone +49 36628 66 33 00

Fax +49 36628 66 33 55

E-mail [info@uniprox.de](mailto:info@uniprox.de)

Customer-Service +49 36628 66 33 70

Fax +49 36628 66 33 77

A company of the Bauerfeind Group \_ [www.uniprox.de](http://www.uniprox.de)

Layout and realization by \_ [WWW.TRENDESETZER.EU](http://WWW.TRENDESETZER.EU)

Rev. 3-2018-4\_49000010000018