

# VACUUM INTEGRATED PUMP





## VIP



### DESCRIPTION

The ALPS™ VIP is an air expulsion integrated pump for elevated vacuum for prosthetic socket. The use of VIP provides patients continuous total contact when donning the prosthesis. This is essential for an even weight distribution and to distribute the suspension of the prosthesis over the whole surface of the residual limb. Elevated vacuum created by the VIP stabilizes the residual limb, overcoming all issues connected with the residual limb changing volume during the day.

### ADDED BENEFITS

-  Efficiently deals with stump volume fluctuation which in turn leads to improved prosthesis control and comfort
-  Reduces the intensity of peaks in pressure on the patient's limbs, improving comfort and providing greater freedom of movement
-  Helps preserve residual limb skin health by improving circulation and consequently increasing limb oxygenation
-  Does not affect the cosmetic appearance of the prosthesis

Easy to  
install

Improves  
blood flow


Improves  
comfort

Improves  
walking

Quiet

Light


# AVAILABLE IN 2 DIFFERENT VACUUM LEVELS

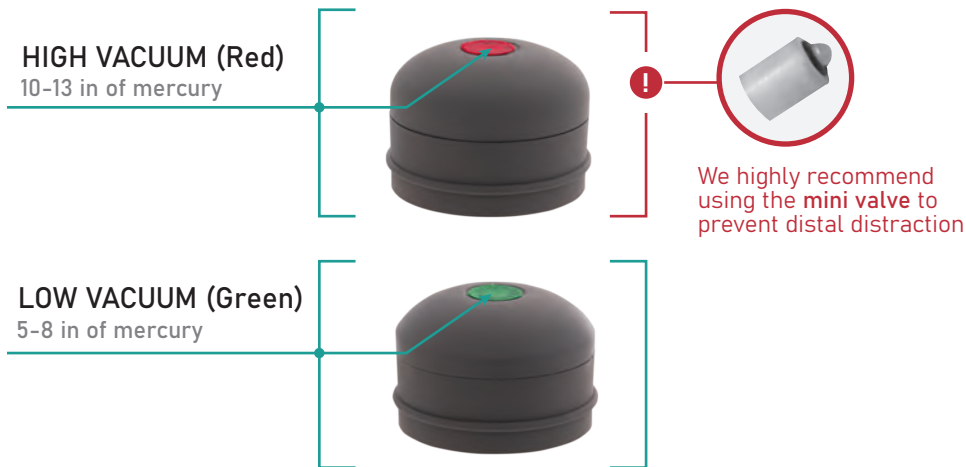


**CHOOSE LEVEL BASED ON:**

- ACTIVITY LEVEL**
- SKIN CONDITION**

## WITHOUT MOUNTING PLATE

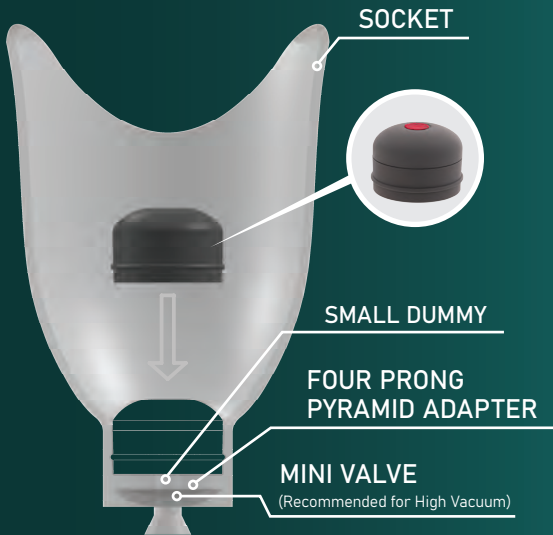

## WITH MOUNTING PLATE

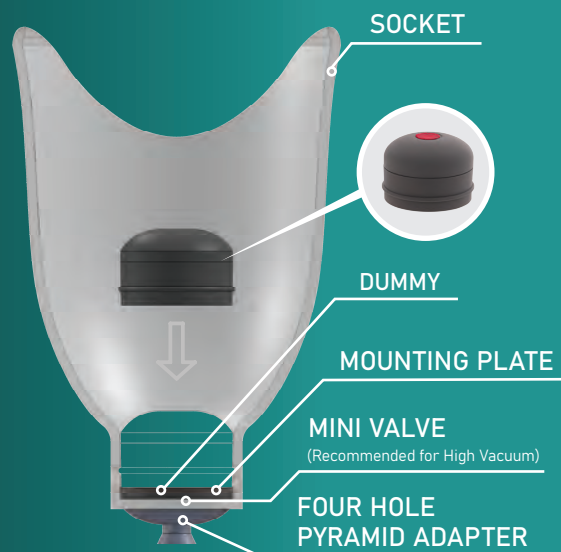

### PRODUCT SPECS

|          | HIGH VACUUM (RED) | LOW VACUUM (GREEN) |
|----------|-------------------|--------------------|
| PART #   | VIP - DH          | VIP - DL           |
| DIAMETER | 66MM              | 66MM               |
| HEIGHT   | 48MM              | 48MM               |
| WEIGHT   | 75G               | 64G                |
| MATERIAL | SILICONE          | SILICONE           |

### USE CODES BELOW WHEN ORDERING

| PART NAME                       | PART #      |
|---------------------------------|-------------|
| DUMMY                           | VE319 - D   |
| SMALL DUMMY                     | VE319 - SD  |
| MOUNTING PLATE                  | VE819 - P   |
| LOCK (SMALL)                    | VIP-L1120-S |
| MINI VALVE (Highly Recommended) | MV-009      |
| MINI VALVE KIT                  | LMV-009     |



T: 727.528.8566 | 800.574.5426 | F: 727.528.8862

www.easyliner.com | info@easyliner.com

2895 42nd Avenue North St. Petersburg, Florida 33714

